


THE MOST STABLE CUSHION FOR WHEELCHAIR USERS


Ischial well - Two compartments
for superior pressure reduction in the IT-coccyx area & extra positioning possibilities

Trochanter support
firmer density to provide lateral stability | redistributes pressure to help create an ischial well

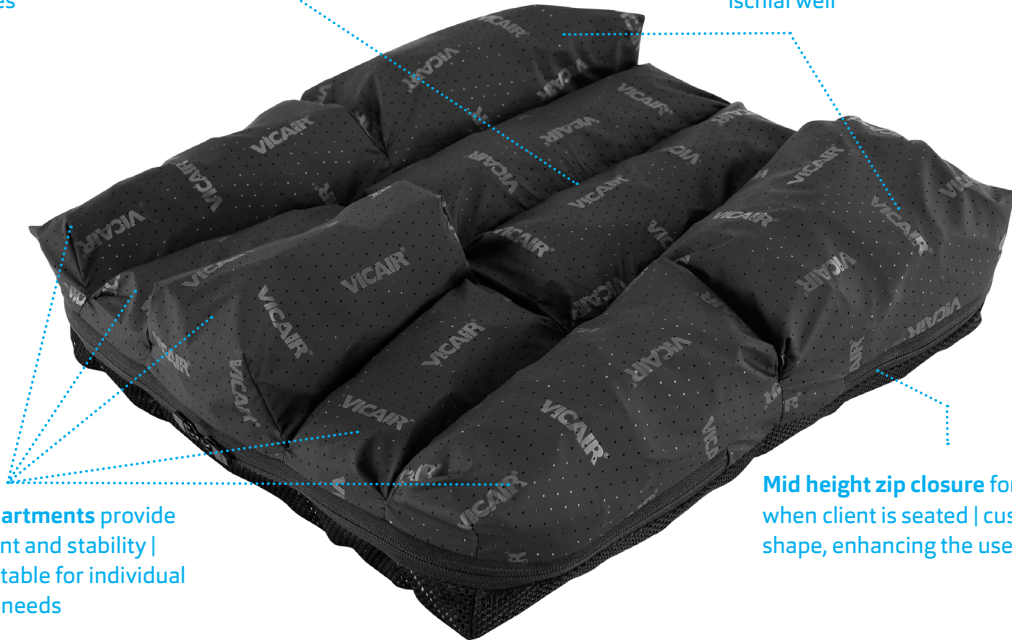

Front compartments provide leg alignment and stability | easily adjustable for individual positioning needs

Mid height zip closure for easy adjustment when client is seated | cushion keeps its shape, enhancing the user's skin protection

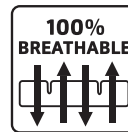

Nine adjustable compartments for optimal stability

Expanding the user's reach

Extra positioning of the thighs

The Vicair Vector O2 cushion offers enhanced stability to the wheelchair user. This results in improved functionality of the upper body and expands the reach of the user's arms. The cushion's unique compartment layout reduces pressure in the vulnerable area around the sitting bones and tailbone. The Vector O2 provides excellent leg guidance, optimal skin protection and extra pelvic stability.

The cushion is 100% machine washable and breathable, for optimal hygiene and perfect temperature and moisture regulation.

PRODUCT PERFORMANCE

Skin protection / Anti decubitus	▲▲▲▲▲
Pressure redistribution	▲▲▲▲▲
Shear force reduction	▲▲▲▲△
Microclimate control	▲▲▲▲▲
Impact & vibration damping	▲▲▲▲▲

Positioning

Lateral stability	▲▲▲▲▲
Frontal stability	▲▲▲▲▲
Sliding reduction	▲▲▲▲▲
Allows manual adjustments	▲▲▲▲▲
Allows easy transfers	▲▲▲▲△

Maintenance & Safety

Durability	▲▲▲▲▲
Reliability	▲▲▲▲▲
Functionality when misused	▲▲▲▲△
Low maintenance level	▲▲▲▲▲
Level of hygiene	▲▲▲▲▲

Functionality & Application

Comfort	▲▲▲▲▲
Low cushion weight	▲▲▲▲▲
Available size range	▲▲▲▲▲
Maximum user weight	▲▲▲▲▲
Low seat to floor height	▲▲▲▲▲
Auto adjustments	▲▲▲▲△

USED WITH

Spinal cord injuries, neuromuscular diseases, hypermobility syndrome, Ehlers Danlos Syndrome (EDS)

TECHNICAL SPECIFICATIONS

Height	6 and 10 cm
Average weight	725 grams
Max. user weight	250 kg
Machine washable	60°C max.

CUSHION SIZES

■ Standard (Sizes >53 cm height 10 only)
■ On request

