



COMFORT AND SKIN PROTECTION IN A TILT-IN-SPACE WHEELCHAIR

Longer compartments to fit into tilt-in-space chairs | follow the body's contours for equal pressure distribution | support IT-coccyx area | managing pelvic asymmetry

Pre-ischial bar limits sliding forward | converts shear into normal forces



Thigh support

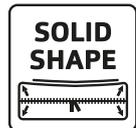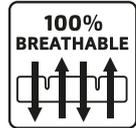

- ▲ Ideal for use in tilt-in-space chair with adjustable seat depth
- ▲ Vicair Adjuster O2 design with deeper rear end compartments
- ▲ Optimal adjustment to the body's contours for equal pressure distribution

The Vicair Multifunctional O2 consists of 5 individual compartments, designed to contour under loading. The cushion features a special 8 cm thick design with elongated rear end compartments for optimal use in tilt-in-space wheelchairs. 8 cm thickness is thick enough to protect the skin and low enough for optimal accessibility under tables, counter tops, desks etc. Each compartment is designed to influence a specific area of the pelvis, striking a unique balance between preserving skin and underlying tissue at the ischials, coccyx and trochanters, positioning the pelvis and supporting the femurs. For specific needs the filling of each of the 5 compartments can be easily adjusted, allowing you to change the density of each compartment.

**AVAILABLE COVERS**  
Vicair O2 Top Cover, Comfair Cover and Incotec Cover.

**PRODUCT PERFORMANCE**

**Skin protection / Anti decubitus**

Pressure redistribution	▲▲▲▲▲
Shear force reduction	▲▲▲▲▲
Microclimate control	▲▲▲▲▲
Impact & vibration damping	▲▲▲▲▲

**Positioning**

Lateral stability	▲▲▲▲△
Frontal stability	▲▲▲▲▲
Sliding reduction	▲▲▲▲▲
Allows manual adjustments	▲▲▲▲▲
Increasing endurance	▲▲▲▲△

**Maintenance & Safety**

Durability	▲▲▲▲▲
Reliability	▲▲▲▲▲
Functionality when misused	▲▲▲▲△
Easy to maintain	▲▲▲▲▲
Level of hygiene	▲▲▲▲▲

**Functionality & Application**

Comfort	▲▲▲▲▲
Ultra light weight	▲▲▲▲▲
Allows easy transfers	▲▲▲△△
Available standard sizes	▲▲△△△
Maximum user weight	▲▲▲▲△

**TECHNICAL SPECIFICATIONS**

Thickness	8 cm
Average weight	850 grams
Max. user weight	250 kg
Machine washable	60°C max.

**CUSHION SIZES**

cm	inches
37x52	14¾x20¾
40x52	16x20¾
45x52	18x20¾
50x52	20x20¾

Sizes in width x depth