

## OPTIMAL CUSHION FOR POSITIONING

### Follows the body's contours

for equal pressure distribution | support IT-coccyx area | managing pelvic asymmetry

### Pre-ischial bar limits

sliding forward | converts shear into normal forces

Mid height zip closure for easy adjustment when clients are seated | cushion keeps its shape, enhancing the user's skin protection

Thigh support



**PRESSURE INJURY TREATMENT PLAN**

**100% BREATHABLE**

**SOLID SHAPE**

**MACHINE WASHABLE**

**READY FOR USE**

- ▲ Minimizes sliding and shear forces
- ▲ Individual solution for asymmetry
- ▲ Optimal adjustment to the body's contours for equal pressure distribution

The Vicair Adjuster 12 02 consists of 6 individual compartments, designed to contour under loading. It automatically manages obliquities and is individually adjustable to support more extreme asymmetry. The thickness of the cushion ensures perfect positioning. Each compartment is designed to influence a specific area of the pelvis, striking a unique balance between preserving skin and underlying tissue at the ischials, coccyx and trochanters, positioning the pelvis and supporting the femurs. For specific needs the filling of each of the 6 compartments can be easily adjusted, allowing you to change the density of each compartment.

### PRODUCT PERFORMANCE

#### Skin protection / Anti decubitus

Pressure redistribution	▲▲▲▲▲
Shear force reduction	▲▲▲▲▲
Microclimate control	▲▲▲▲▲
Impact & vibration damping	▲▲▲▲▲

#### Positioning

Lateral stability	▲▲▲▲△
Frontal stability	▲▲▲▲▲
Sliding reduction	▲▲▲▲▲
Allows manual adjustments	▲▲▲▲▲
Increasing endurance	▲▲▲▲△

#### Maintenance & Safety

Durability	▲▲▲▲▲
Reliability	▲▲▲▲▲
Functionality when misused	▲▲▲▲△
Easy to maintain	▲▲▲▲▲
Level of hygiene	▲▲▲▲▲

#### Functionality & Application


Comfort	▲▲▲▲▲
Ultra lightweight	▲▲▲▲▲
Allows easy transfers	▲▲△△△
Available standard sizes	▲▲▲▲▲
Maximum user weight	▲▲▲▲△

### TECHNICAL SPECIFICATIONS

Thickness	12 cm
Average weight	900 grams
Max. user weight	250 kg
Machine washable	60°C max.

### CUSHION SIZES

- Standard
- Special Size
- Custom Size

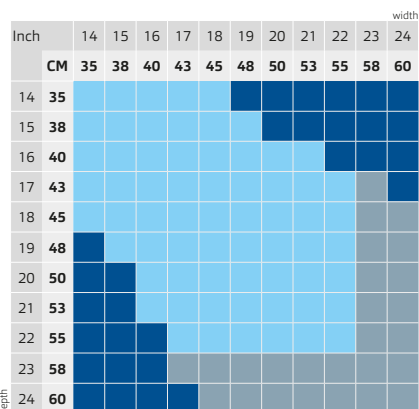

### AVAILABLE COVERS

Vicair O2 Top Cover, Comfair Cover, Incotec Cover and Reversible Cover.