
## VICAIR<sup>®</sup> JUNIOR VECTOR O2

## THE OPTIMAL PEDIATRIC CUSHION

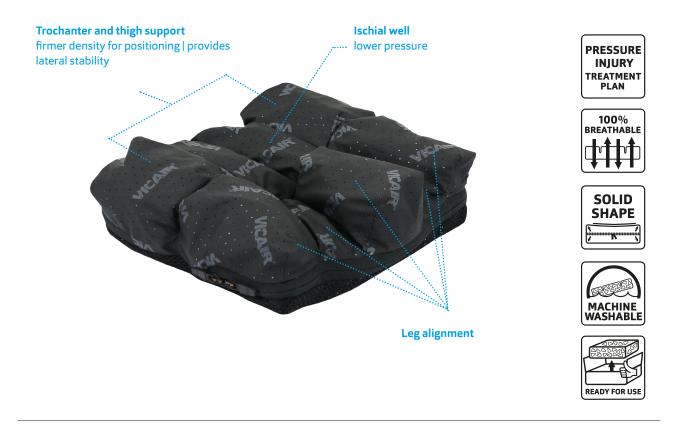

- Especially designed for children
  - Filled with small-sized SmartCells for optimal pressure redistribution
- Built to provide carefree mobility

The Vicair Junior Vector O2 is a precontoured cushion which consists of 9 individual compartments. This cushion is designed especially for smaller people and pediatric use. Unique are its special, smaller sized SmartCells that safeguard pressure management and positioning properties even in the smallest of cushion sizes. Due to the anthropometrics of the pelvis of the target group, the thickness of only 6 cm is excellent for individual positioning and skin protection. Each compartment is designed to influence a specific area of the pelvis, striking a unique balance between preserving skin at the ischials and coccyx by redistributing pressures, stabilizing the pelvis and positioning of the femurs. For specific needs the 9 compartments can be easily adjusted.

## AVAILABLE COVERS

Vicair O2 Top Cover, Comfair Cover and Incotec Cover.

PRODUCT PERFORMANCE	
Skin protection / Anti decubitus	5
Pressure redistribution	
Shear force reduction	
Microclimate control	
Impact & vibration damping	
Positioning	
Lateral stability	
Frontal stability	
Sliding reduction	
Allows manual adjustments	
Increasing endurance	
Maintenance & Safety	
Durability	
Reliability	
Functionality when misused	
Easy to maintain	
Level of hygiene	
Functionality & Application	
Comfort	
Ultra lightweight	
Allows easy transfers	
Available standard sizes	
Maximum user weight	

TECHNICAL SPECIFICATIONS					
Thickness	6 cm				
Average weight	450 grams				
Max. user weight	100 kg				
Machine washable	60°C max.				

CU	CUSHION SIZES							
	Sta	nda	ırd					
On request								
						width		
Inch		10	11	12	13	14		
	CM	25	20	20		25		

		СМ	25	28	30	33	35
	10	25					
	11	28					
	12	30					
	13	33					
Indar	14	35					

070069-0JVR-EN-REV01/2024